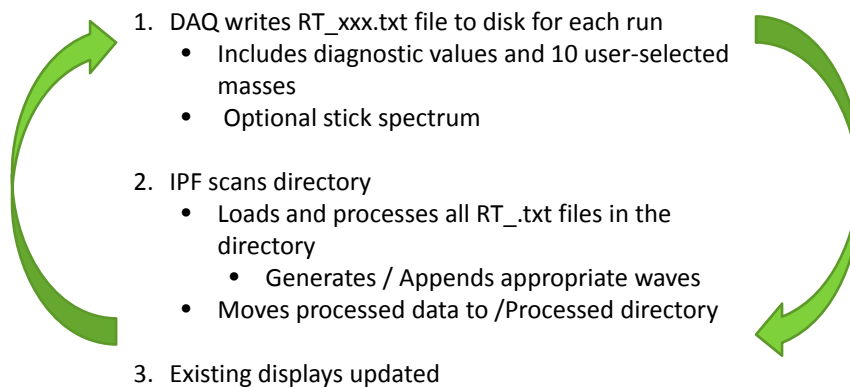


Realtime .txt (aka Quickview)


Joel Kimmel
jkimmel@aerodyne.com

AMS Users Meeting
Toronto
31-Oct-2009

Realtime.txt Process

1. DAQ writes RT_xxx.txt file to disk for each run
 - Includes diagnostic values and 10 user-selected masses
 - Optional stick spectrum
 2. IPF scans directory
 - Loads and processes all RT_.txt files in the directory
 - Generates / Appends appropriate waves
 - Moves processed data to /Processed directory
 3. Existing displays updated
- 

Enables customized real time processing and displays

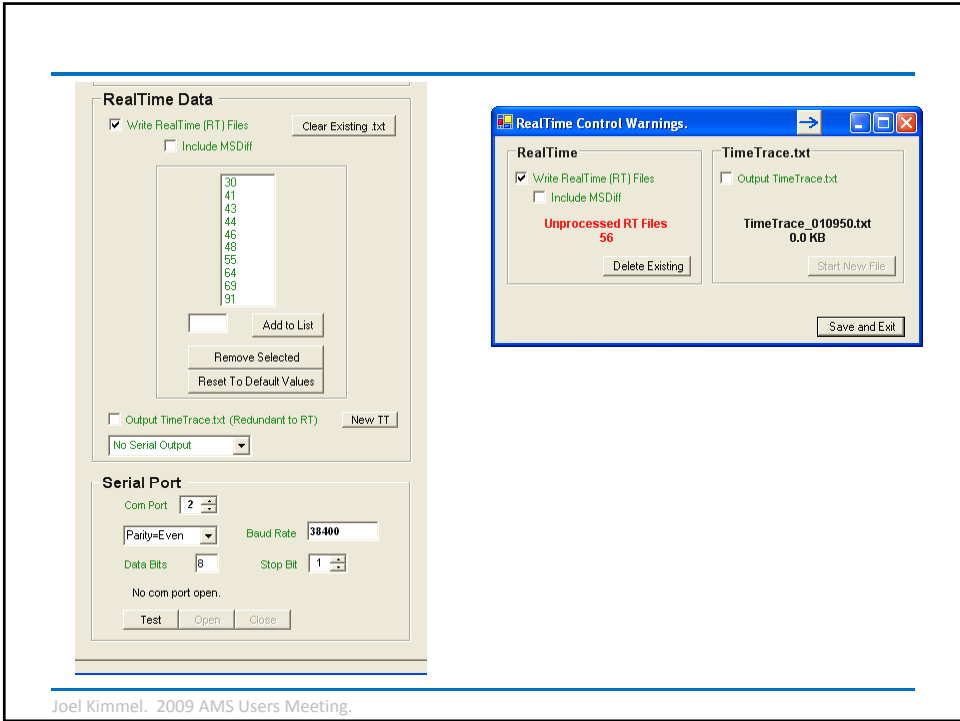

Joel Kimmel. 2009 AMS Users Meeting.